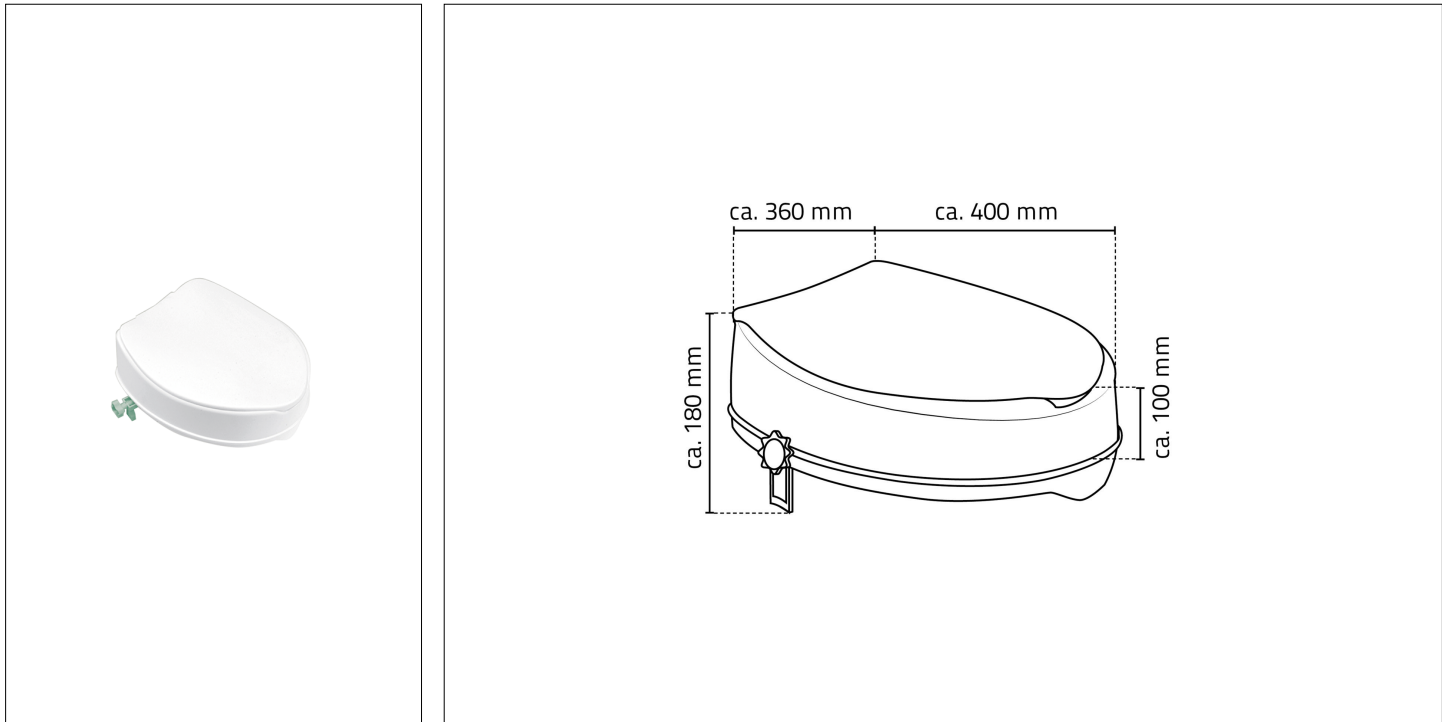




Product code	515	Product Family	SEDILI WC	Collection	Sanitaryware accessories
Description	Raised toilet seat Sedile universale H= 100 mm con apertura frontale e coperchio. Fixing via side clamps. Raised, H=100 mm				
Certifications	TUV UNI EN ISO 9001 system certification.				CE
Prestazioni	Load capacity 150 kg.				
Features	H100 universal seat with cover and front opening, in low toxicity polypropylene. For people who have difficulty sitting or getting up when using the toilet, such as the elderly and post-operative users. The seating area is ergonomic and reduces pressure offering greater comfort. The slight forward inclination of the seat helps maintain correct posture and reduces flexion of the hip and knee joint. Seat characterized by two openings, one in the front and one in the back which allow better access for personal hygiene. Adatto alla maggior parte dei water. Portata massima 150 kg. Dimensioni (indicative): Lunghezza 40 cm Larghezza 36 cm. Superficie lucida facile da pulire e resistente alle macchie ed agli odori. La presenza del coperchio aiuta a mantenere l'igiene del prodotto. L'alzawater si fissa sul water in maniera stabile e sicura tramite un meccanismo di bloccaggio composto da due staffe in nylon laterali regolabili in larghezza. Il meccanismo di fissaggio offre resistenza e affidabilità; il design ergonomico della manopola permette una facile e veloce applicazione e rimozione dell'alzawater anche a coloro che abbiano disfunzioni alle articolazioni della mano.				
Finishes	White colour.				
Dimensioni	410 x 175 x 370				
Installation	Secure to the vase using the appropriate clamps. The manufacturer is not responsible for the installation of the product. Periodically check that the fixing is solid and stable. Keep the product properly clean.				

Maintenance	Periodically check the correct fastening of the accessories. Carry out adequate cleaning with soap and water, do not use aggressive, abrasive, corrosive products, alcohol, acids, acetone, thinner, etc. Don't scratch. Disinfect periodically. Do not use the product for purposes other than those intended. After cleaning, leave the lid and seat raised to prevent them from being damaged by the vapors of the sanitaryware cleaning products.
Packaging	Imballo in scatola di cartone 410 x 175 x 370 mm / Imballo Peso: 2 kg
Weight	2 kg